
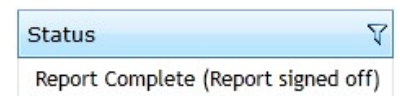


OVERVIEW

The recovery screen within MEDILOGIK EMS™ allows staff to capture information about the patient’s episode of care once the report has been signed off by the Endoscopist.

ACCESS THE PATIENT RECORD



1. Log in to EMS.
2. Click the “Daily List” panel. If this is not visible, select **Home / Panels / Daily List** to activate the panel.
3. Ensure that the patients report has been signed off. This does not affect the ability to complete the recovery screen but it is important for record keeping.
4. Click the name of the patient to be taken to the latest step in their care pathway.

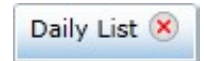


RECOVERY

NOTE: An image of the recovery screen is shown overleaf.

1. Record the recovery start time.
2. In the “Recovery” screen you can record information about any post procedure medication or adverse events.
Click [**Nursing Observations**] to open the observations window. Click [**New Observation**] to enter the details of any observations. Click [**OK**] to complete.
3. Repeat the addition of observations as necessary.
4. Record the recovery finish time.
5. Click [**Discharge**] to progress.

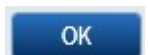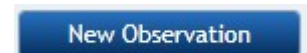


THE RECOVERY SCREEN

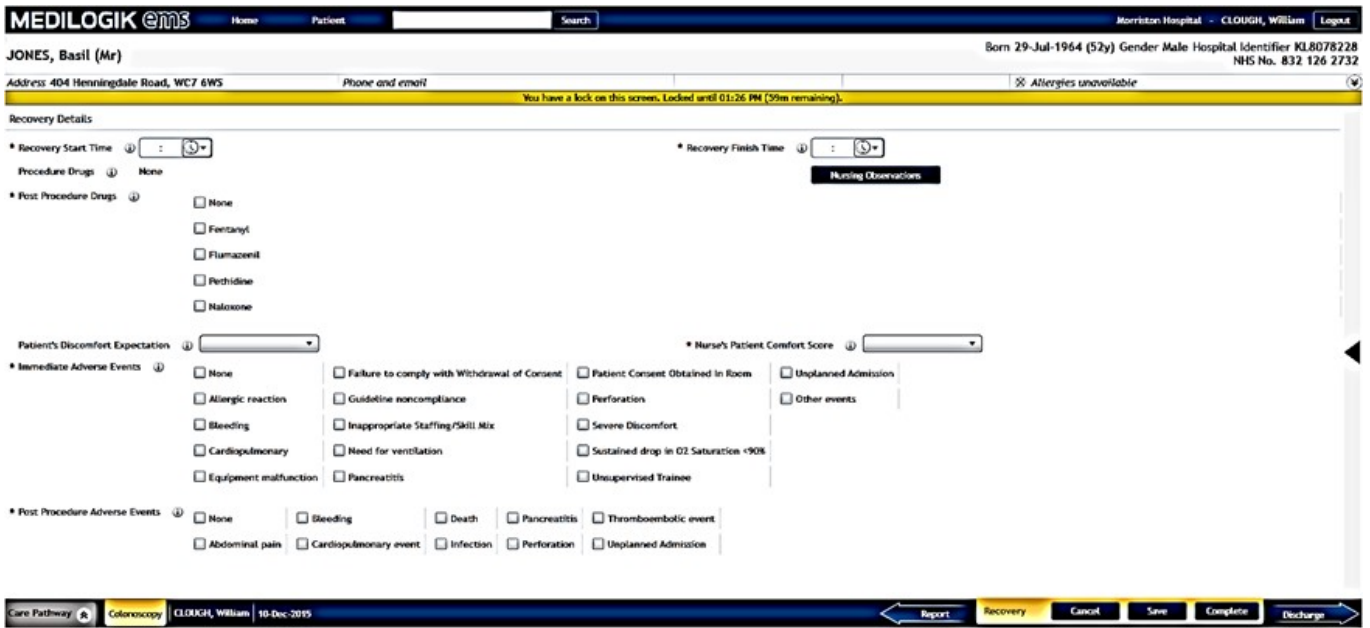

DISCHARGING A PATIENT

Following the period of recovery the patient will be available to be discharged.

1. Record the patient's discharge time.
2. Confirm if any further information has been given to the patient.
3. If necessary, users can record comments relating to the discharge.
4. Click **[Complete]** to complete the episode of care.

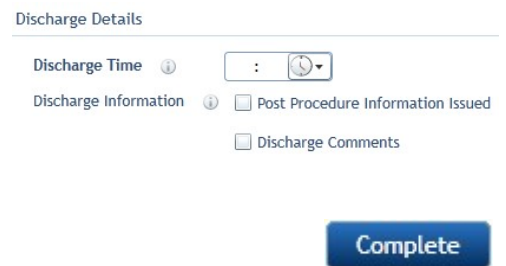

STATUS

1. Having discharged the patient, click **[Home]** to return to the "Daily List" panel.
2. The screen should refresh automatically and the patient's episode status should now show "Discharge Complete (Signed Off)". This confirms the completion of the entry.

